100% DOUBLE QUARTZ GLASS ELECTRODE



100% DOUBLE QUARTZ GLASS ELECTRODE

This Medical Ozone Generator uses DOUBLE PURE QUARTZ WALLS between the electrodes, a new patent design that offers Ultrapure Ozone and constant stable ozone concentration output. (see Page 6 for patent model)



This ozone charmber is only with pure quartz, has no any metal nor plastic, rubber, etc. at the ozone generating area. It makes truly Contaminant-Free ozone for treatments.

Key features:


- Ultra purity of ozone for medical use
- ➤ High accuracy of ozone output
- Constant Stable ozone concentration
- Multi-levels for options
- ➤ Built-in Ozone Destructor
- with Patents



100% DOUBLE QUARTZ GLASS ELECTRODE

Model No.	MOP0.3-AD
Ozone Concentrations:	$2 \sim 108 \text{ ug/mL}$ ($1 \text{ ug/mL} = 1 \text{ gamma} = 1 \text{ mg/L}$) Adjustable
Connection of Oxygen Inlet, Ozone Outlet & Ozone Destructor Inlet	Standard Luer Lock Connector
Cooling System	4000 RPM, turbo brushless cooling fan
Ozone Module	100% Double Quartz Glass Electrode (Patent Design)
Life	> 10,000 hours
Input Voltage	12V DC (100-240V AC to 12V DC adaptor)
Rated Power	< 12W
Operating Temperature	-10 °C ~ 40 °C
Storage Temperature	-20°C ~ 70°C
Product Size(L x W x H)	228*176*80 mm
Package Box Size(L x W x H)	330 *230*150 mm
Gross weight	1.5 kg (if with oxygen tank regulator: 1.9 kg)

100% DOUBLE QUARTZ GLASS ELECTRODE

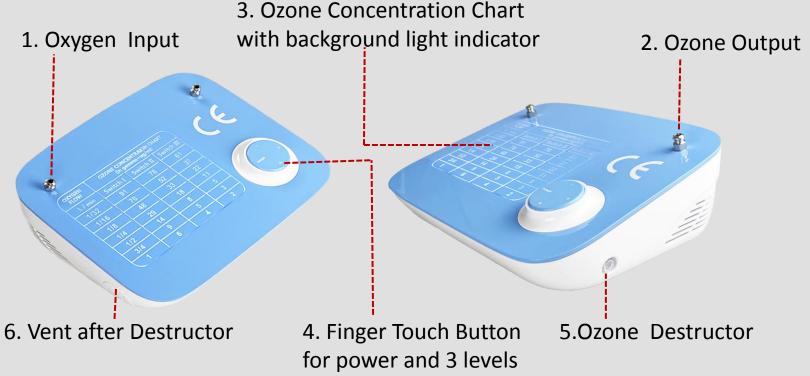

Accessories:



1pc *
Power Adaptor



1pc *
3 way connector



3pcs * male luer connector /



1.5 meter * silicon tubing

100% DOUBLE QUARTZ GLASS ELECTRODE








L/min		OZONE CONCENTRAION CHART (in gamma=ug/ml) Switch I Switch II Switch III			2
1/32	91	78	61		
1/16	70	52	37		
1/8	48	33	22		
1/4	29	18	11		
1/2	14	8	5	**************************************	
3/4	9	5	3		
1	6	4	2		



100% DOUBLE QUARTZ GLASS ELECTRODE

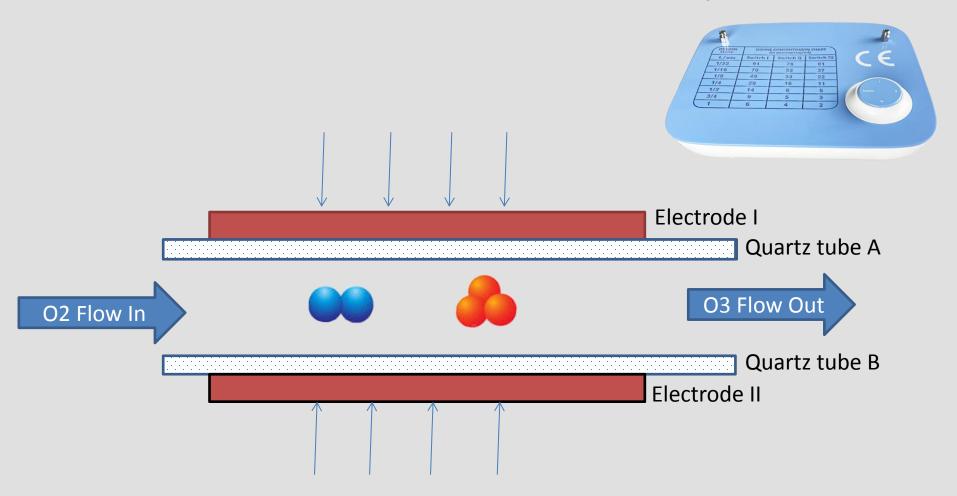

The material of ozone generating area: ONLY PURE QUARTZ GLASS contacted the gas flow

100% DOUBLE QUARTZ GLASS ELECTRODE

100% DOUBLE QUARTZ GLASS ELECTRODE

Made for medical use

Ozone Concentration Charts

Electrode by 100% Pure Quartz Glass

Double Quartz Walls with patents

Ideal for different ozone treatments

C Compliance with CE ROHS

Available for global travelling

Latest technology

)				6
XYGEN FLOW	OZONE (CONCENTRAIO gamma=ug/r		
L/min	Switch I	Switch II	Switch III	
1/32	108	94	76	66
1/16	81	61	44	
1/8	51	35	24	
1/4	29	19	12	
1/2	15	9	6	
3/4	10	6	3	
1	7	4	2	

Notice: The "Ozone Concetnration Chart" might be slightly different from the above picture.

Each machine will be tested and calibrated before shipment







Ozone joint injection



Ozone blood treatment



Ozone Rectal Insufflation



Gynecological ozone therapy



Ozone for dental



Ozone Drinking Water



Ozone Ear Insufflation



Ozone for Wound Healing



Ozone for making Ozone Oil



Ozone Sauna



Ozone for Lyme disease treatment



Ozone for Cancer treatment

Recommended applicatoins

Product Package







Standard Package includes:

- Ozone device * 1pc
- Power adaptor * 1pc
- Silicon tubing * 1.5 meters
- 3 way connector * 1pc
- Luer connector * 3pcs

Relative Parts for ozone treatments



O2 regulator CGA870



O2 regulator CGA540



O2 regulator bullnose



O2 regulator Brazil



O2 regulator Germany



O3 destructor



power adaptor



Glass cup



SS#316 bubbler



Silicon tubing



flask 500mL



flask 1000mL


Relative Parts for ozone treatments










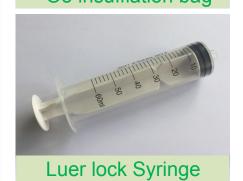

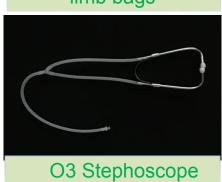













Silicon tubing